

Providing you with gaseous chlorine dioxide solutions for your decontamination and sterilization needs



Distribution Plate (Part #: CSI-ROOM015)

Description

The Distribution Plate provides a method to connect ClorDiSys Solutions, Inc. gas generators to multiple rooms from one central location for quick easy connections. It is mounted at a central location and tubes and cabling run from this panel to the ceiling plate at each individual room to be decontaminated. It is a custom part and can be made for 1-20 different rooms.

The plate allows for quick valved tubing connections and quick connection of electrical connectors to measure and control the RH to either the Mindox or Cloridox systems

This part is used in conjunction with CSI-ROOM014 which is mounted in the room.

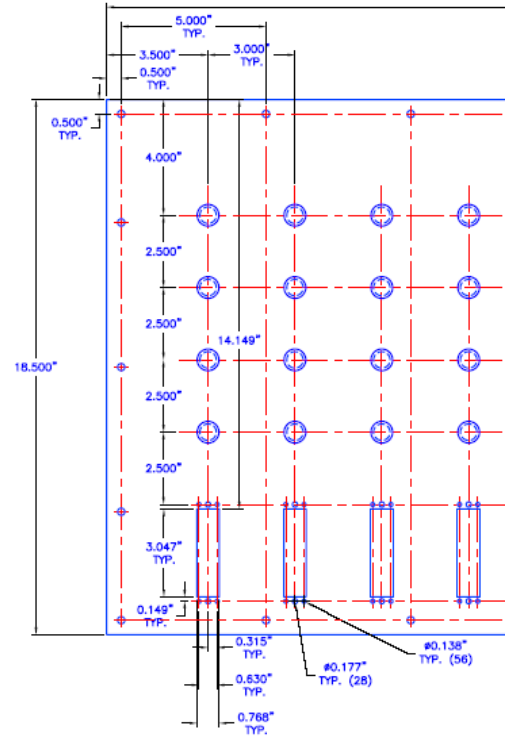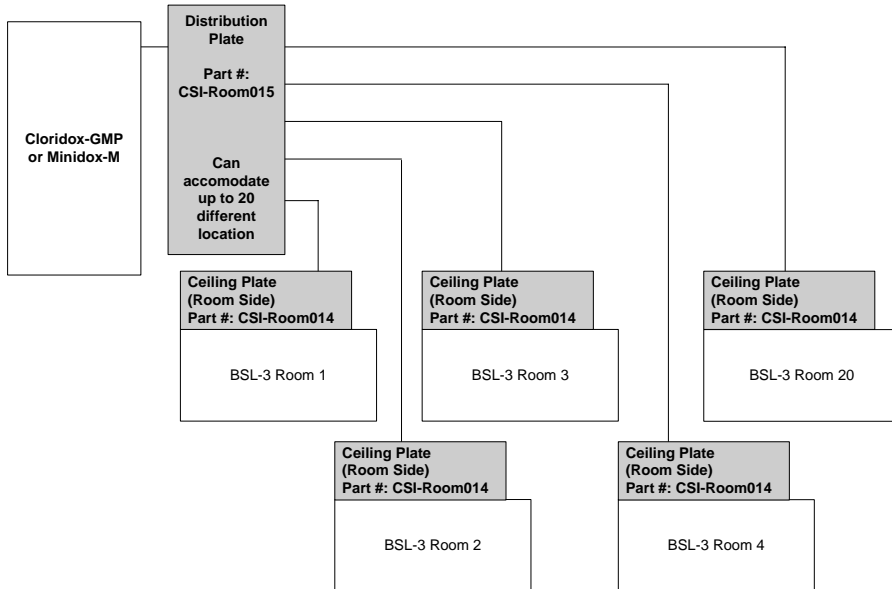

Distribution Plate (Part #: CSI-ROOM015)
Unit Dimensions for 14 chambers: 16.0" x 56.0"
Unit Weight: 15 lbs

